

Genotyping using Extract-Amp RED PCR Kit



Genotyping is a technology used for detection of minor genetic differences. Genotyping determines differences in genetic material by comparing one DNA sequence to another DNA sequence or a reference. Genotyping is used in research and diagnostic medicine.

Extract-Amp RED PCR Kit from Apex is the ideal, fast and convenient solution for genotyping. The performance of Extract-Amp RED PCR Kit has been tested using different mammalian matrices and also against similar genotyping PCR kits from other suppliers. Furthermore, the performance and handling time of the included Apex DNA Extraction Solution was compared to similar fast DNA extraction solutions from other suppliers.

Features

- Genotyping using mammalian tissue e.g. mouse ear or tail*
- PCR-ready DNA in 8 minutes
- Minimal handling
- Reliable PCR results
- Red dye for direct gel loading and visualisation of pipetting
- Scalable set-up
- Automation-friendly
- Results obtained 1-2 hours after sample collection

*Other samples such as saliva, bacteria, plants, bird feathers, fish fins and various mammalian tissues can be used.

DNA is extracted in just 8 minutes from many different mammalian sample types e.g. mouse tails, ears and saliva. Depending on the sample size, the DNA extraction is performed in either PCR tubes or 1.5 ml tubes using a thermocycler or heating block, respectively. The one-reagent one-tube extraction set-up is easily scaled for e.g. automation on robotic platforms.

Extract-Amp RED PCR Kit combines the user-friendly DNA extraction of Apex DNA Extraction Solution with the convenience and excellent PCR performance of 2.0X Taq RED Master Mix. The non-toxic Apex DNA Extraction Solution is designed for rapid and efficient extraction of PCR-ready for genotyping.

Efficient amplification of DNA extracted from different types of mammalian tissues

Extract-Amp RED PCR Kit was used to extract and amplify genomic DNA from five different mouse tissues (kidney, muscle, lung, liver and tail), chicken muscle tissue and human saliva, (fig. 2). For DNA extraction, 0.5 -10 mg tissue or 20 µl saliva was added to 100 µl of Apex DNA Extraction Solution. The extracted DNA was then amplified using the 2.0X Taq RED Master Mix and species specific primers. The results show that the Extract-Amp RED PCR Kit can be used for many different mammalian samples and also saliva.



Figure 2. Amplification of DNA from various samples. Extract-Amp RED PCR Kit was used to extract and amplify genomic DNA from various mammalian tissues. M: DNA marker ECON Low DNA Ladder. Lane 1-5: Different mouse tissues as depicted, GAPDH (266 bp). Lane 6: Chicken muscle tissue, HRPT1 (245 bp) and Lane 7: Human saliva, DMD17 (415 bp).

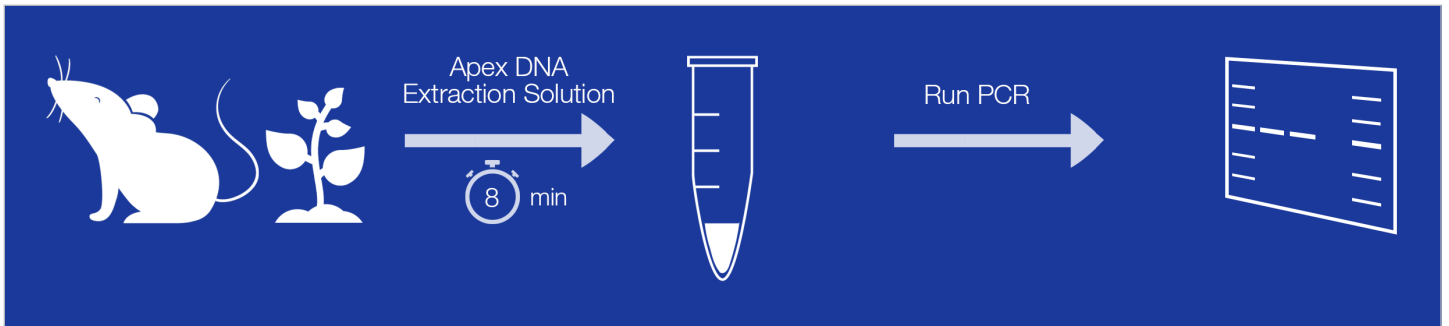


Figure 1. Illustration of the genotyping workflow, when using Extract-Amp RED PCR Kit. The PCR-ready DNA is extracted in 8 minutes (6 min at 65 °C + 2 min at 98 °C), from mouse tail tissue or other mammalian tissues, using the Apex DNA Extraction Solution. The extracted DNA is ready for PCR without further handling such as vortex, centrifugation or dilutions. The extracted DNA is amplified using 2.0X Taq RED Master Mix.

Extract-Amp RED PCR Kit

Flexible sample handling

The flexibility of sample handling was tested by adding varying amounts (mg) of chicken muscle tissues to 100 µl of Apex DNA Extraction Solution. The extracted DNA was subsequently amplified using 2.0X Taq RED Master Mix. The results show that the Apex DNA Extraction Solution protocol provides the user with a high degree of flexibility over a wide range of applied sample amounts, (fig. 3).

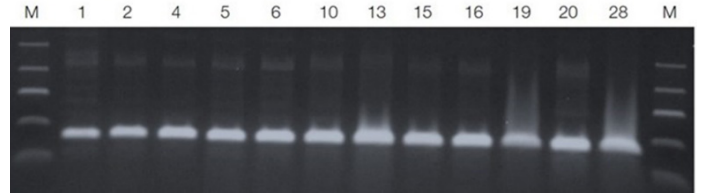


Figure 3. Flexibility of Apex DNA Extraction Solution protocol. M: DNA marker ECON Low DNA Ladder. Lanes 1-28: Varying amounts (mg) of chicken muscle tissue, HRPT1 (245 bp) as specified on the top of each lane.

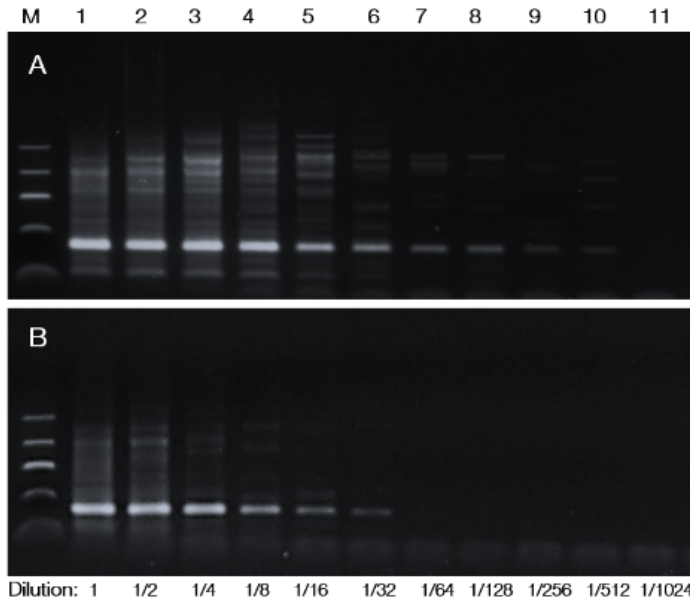


Figure 4. The performance of Extract-Amp RED PCR Kit when using either the standard 3-step PCR protocol, panel A, or the Fast 3-step PCR protocol, panel B. Extract-Amp RED PCR Kit was used to extract and amplify genomic DNA from mouse tail tissue. M: DNA marker ECON Low DNA Ladder. Lane 1-11 Amplification of a two-fold serial dilution of the extracted mouse genomic DNA, mouse GADPH (266 bp).

Fast genotyping results

The Extract-Amp RED PCR Kit was used to extract and amplify 2-fold dilutions of DNA from mouse tail tissue. The amplification was performed using either a standard 3-step PCR protocol, panel A or a fast 3-step PCR protocol, panel B, (fig. 4). Specific PCR results are obtained using both standard and fast 3-step protocol. We have recently shown that it is possible to save ~1 hour of PCR run time using fast 2-step-protocol and still get specific PCR amplicons of high yields. When using Extract-Amp RED PCR Kit, results (from sample to result) are obtained in 70 to 210 minutes, depending on the combination of applied DNA extraction method as well as applied PCR protocol, (fig. 5).

Storage of DNA extracts

DNA extracts can be stored at -20 °C for one week.
Long term storage at -80 °C.

Standard 3-step PCR protocol:

Step	Temperature	Time	Cycles
Initial	95 °C	5 min	1
Denaturation	95 °C	30 sec	35
Annealing	55 °C	30 sec	
Elongation	72 °C	30 sec	
End	4 °C	∞	1

Fast 3-step PCR protocol:

Step	Temperature	Time	Cycles
Initial	95 °C	3 min	1
Denaturation	95 °C	15 sec	35
Annealing	55 °C	15 sec	
Elongation	72 °C	20 sec	
End	4 °C	∞	1

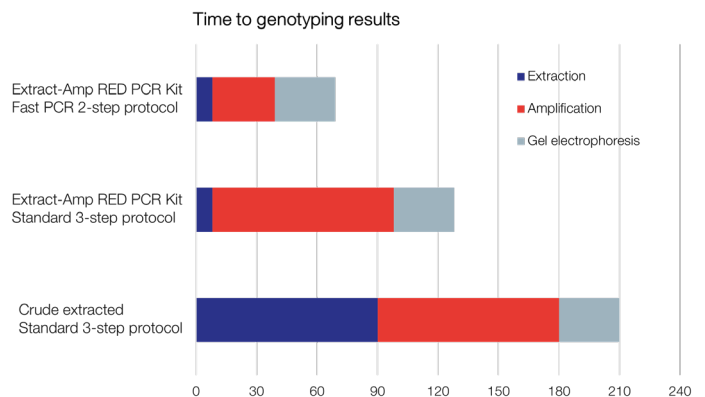


Figure 5. Total time (in minutes) from sample to result required for DNA extraction, amplification and gel electrophoresis when using Extract-Amp RED PCR Kit with either fast 2-step PCR protocol* or standard PCR protocol. Furthermore, time to result of the standard crude extraction protocol** is also included.

*Application note can be downloaded from Apex Taq DNA Polymerase product pages.
**(Truett GE et al. 2000. Biotechniques 29 (1): 52-54.)

Extract-Amp RED PCR Kit

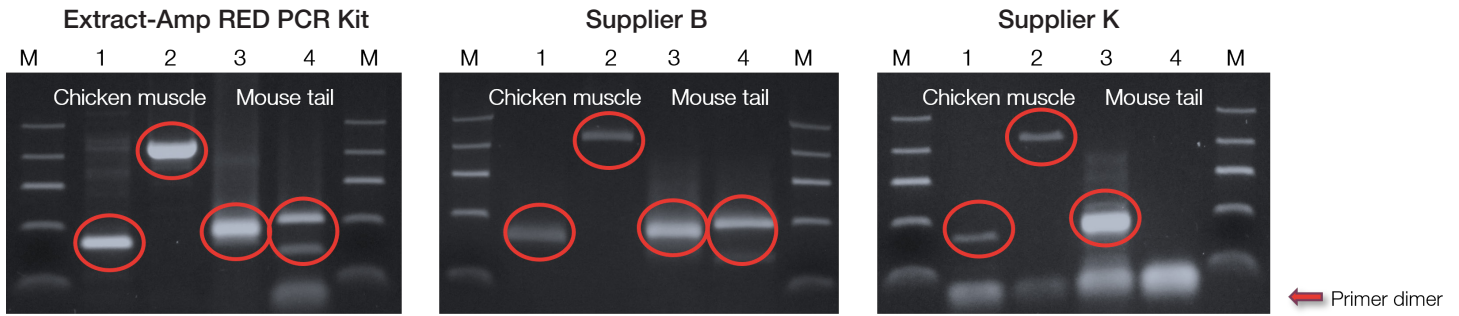


Figure 6. Performance of Extract-Amp RED PCR Kit was compared to similar genotyping PCR kits from supplier B and K. DNA was extracted from chicken muscle (lane 1-2) and mouse tail (lane 3-4) and amplified using DNA extraction and PCR reagents provided in the respective genotyping PCR Kits. M: DNA marker ECON Low DNA Ladder. Lane 1: Chicken HRTP1 (245 bp), Lane 2: Chicken GADPH (775 bp), Lane 3: Mouse GADPH (265 bp) and Lane 4: Mouse B-actin (318 bp). Correct/expected amplicons are encircled.

Performance of Extract-Amp RED PCR Kit compared to two equivalent PCR Kits

DNA was extracted from chicken muscle or mouse tail. Both DNA extracts were tested with two different primer sets, respectively, (fig. 6). Each extraction and amplification were conducted according to the supplier manuals and using the DNA extraction reagents and PCR master mixes provided in the respective kits. The result shows that the Extract-Amp RED PCR Kit performs equally well or better than supplier B and L on the four DNA targets tested. The applied extraction protocol from supplier B and K are very similar and less user-friendly than the Apex DNA Extraction Solution protocol, (fig. 8).

Performance of Apex DNA Extraction Solution compared to an equivalent DNA extraction solution

The performance of Apex DNA Extraction Solution was compared to similar DNA extraction solution from supplier L, (fig. 7). DNA was extracted from 3 mg of chicken muscle tissue using the respective DNA extraction protocols. A four-fold serial dilution of the two DNA extracts was amplified using 2.0X Taq RED Master Mix. The obtained results indicate that the quality and yield of DNA extract using Apex DNA Extraction Solution is slightly better than that of supplier L. The Apex DNA Extraction Solution protocol provided reliable PCR-ready DNA in 9 minutes which is close to the 10 minutes provided by protocol from supplier L, (fig. 8).

Apex DNA Extraction Solution provides the easiest and fastest DNA extraction protocol

Total handling time of the Apex DNA Extraction Solution protocol was estimated in the laboratory and compared to three similar fast DNA extraction protocols performed according to the manuals from the respective manufacturers, L, K and B, (fig. 8). The total handling time for extraction of DNA using Apex DNA Extraction Solution was faster than the handling time for all three suppliers, but very close to supplier L. Furthermore, Apex DNA Extraction Solution only requires 4 handling steps, which is lower than all the three competitors. Numbers of handling steps are indicated on the bars of the respective competitor in fig. 8.

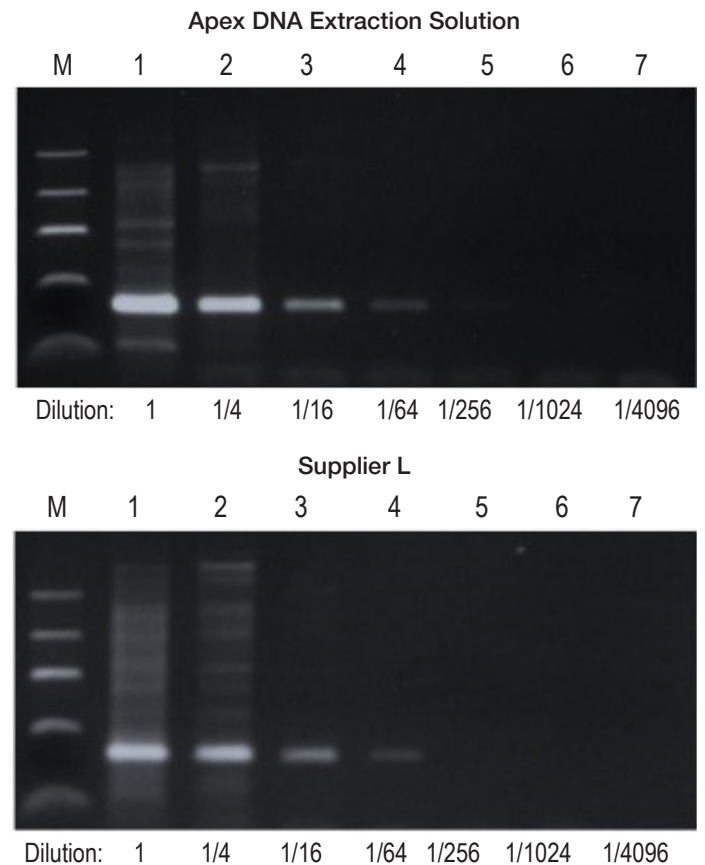


Figure 7. The performance of Apex DNA Extraction Solution was compared to the DNA extraction solution from competitor L. A four-fold serial dilution of DNA extracts was amplified using 2.0X Taq RED Master Mix. M: DNA marker ECON Low DNA Ladder. Lane 1-7: Chicken HRTP1 (245 bp).

Handling time DNA extraction protocols

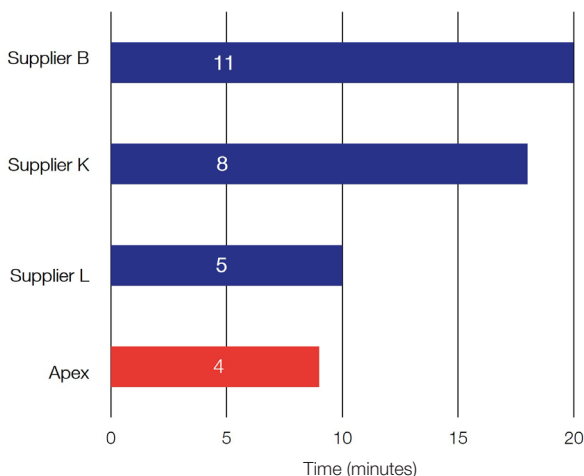


Figure 8. Total handling required for DNA extraction when using either Apex DNA Extraction Solution or similar DNA extraction protocols from supplier L, K and B. Number of handling steps for each protocol is indicated.

Extract-Amp RED PCR Kit

Ordering information

Product	RXN*	Cat #
Apex DNA Extraction Solution	100	42-503
	500	42-503B
<hr/>		
Extract-Amp RED PCR Kit incl. DNA Extraction Solution and 2.0X Taq RED Master Mix	100	42-502
	500	42-502B

* 1 reaction = 100 µl DNA Extraction Solution + 12.5 µl 2.0X Taq RED Master Mix
(final PCR reaction 25 µl)



A Life Science Company

Toll Free: 800.789.5550 Fax: 888.789.0444 Web: www.geneseesci.com Email: support@geneseesci.com