


Stir Bars & Stir Bar Retriever

- Magnetic bar is made from ALNICO, one of the strongest permanent magnets available
- The sintered manufacturing process reduces the chance of decoupling (spin outs), encapsulation in PTFE protects both magnet and what it is stirring
- Chemically inert and easy to clean
- Can withstand temperatures up to 250°C, making them suitable for use with both unheated and heated stirrers
- Available in two surfaces:
 - Low friction octagon shape with a center pivot ring helps create a strong mixing vortex, even with more viscous solutions
 - Smooth shape for better coupling when using concave bottom glassware
- The retrieval bar can recover stir bars even in hot solutions when stirring and mixing is no longer required

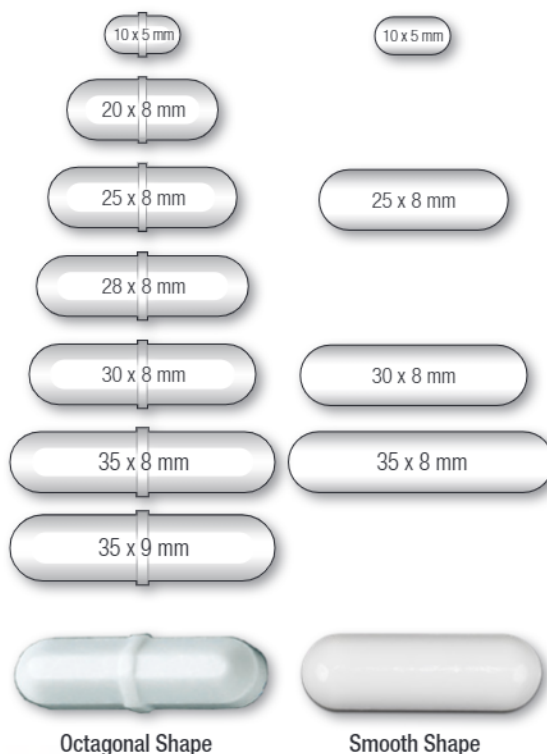

Stir Bar Retriever

Octagonal Shape

Smooth Shape



Stir Bars

Pack of 10

Item No.	Shape	L x W (in)	L x W (mm)
120547	Octagonal with pivot ring	0.4 x 0.2	10 x 5
120548		0.8 x 0.3	20 x 8
120549		1 x 0.3	25 x 8
120550		1.1 x 0.3	28 x 8
120551		1.2 x 0.3	30 x 8
120552		1.4 x 0.4	35 x 9
120894	Smooth	0.4 x 0.2	10 x 5
120895		1 x 0.3	25 x 8
120896		1.2 x 0.3	30 x 8
120897		1.4 x 0.3	35 x 8

Stir Bar Retriever

Pack of 5

Item No.	L (in)	L (cm)
120553	11.8	30