

## TGF $\beta$ 1 inhibitor (TGF $\beta$ 1 binding peptide)

**Code** PG-002

**MW (Acetate):** 2587.90  
**Purity:** >95% (HPLC)

**Amount:** 10  $\mu$ g

**Formulation:** This vial contains 10  $\mu$ g of the titled compound, which has been lyophilized from 0.2  $\mu$ m filtered DMSO solution.


**Reconstitution:** Reconstitute at 50~ 150  $\mu$ M in DMSO  
Centrifuge the vial before opening.

<Example>

100  $\mu$ M solution prepared when 10  $\mu$ g peptide reconstituted by 38.6  $\mu$ L DMSO

**Storage:** Store at -20°C or below for the container vial unopened  
Recommend to use the peptide as quickly as possible after opening and reconstitution.

**Activity:** Measured by its ability to inhibit SBE reporter activity induced by 0.13 nM TGF $\beta$ 1.

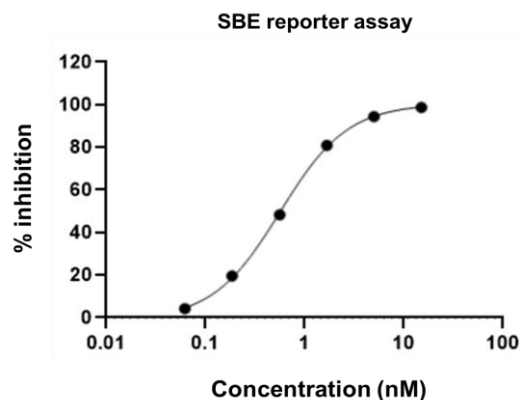

**This product is "For Research Use Only"**